

LUCID`LC

LCD CHART





Key features


- **Google's Android Embeded System** No keyboard! No mouse! Take it with convenient interface only!
- **Multi-control channel provided** Maximum five pieces of the instruments can be operated in the same place.




- **Various visual acuity (VA) chart provided**
Alphabet, numeric, snellen, randolt ring and children charts are provided basically.




- **Various astigmatism chart provided**
Various and useful charts for astigmatism are provided to be applied varied measurement.

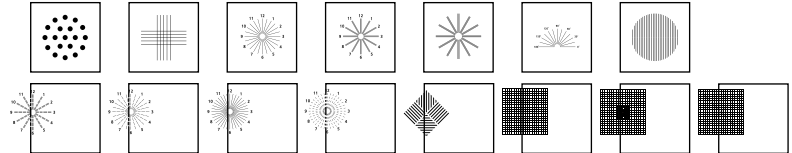

- **Special chart provided**
A lot of special charts for phoria, balance, and stereo are provided as basic options. Any chart provided by operator can be added easily with special tool. For this, consult with your local distributor.

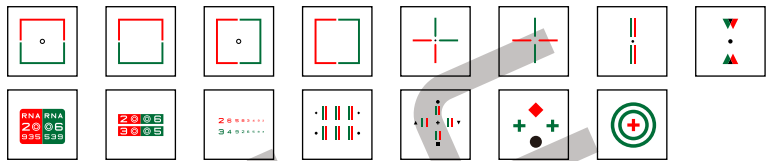

- **Color blindness test**
Basically 23 sheets of high resolution image for color blindness test are provided.

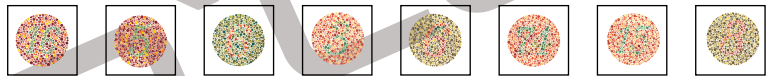

- **ETDRS with contrast sensitivity**
Contrast of ETDRS chart can be adjusted for contrast sensitivity test.



- **Random chart function** Random charts are provided in VA chart in order to prevent patient's memory problem.

- **Reverse function**
It offers to measure the eyesight for post lasek, glaucoma eye or other ophthalmic disease.



Specifications

Display	24 inch WUXGA TFT LCD (1920 x 1200 pixels)
Active display area	(W)518.4 x (H)324.0mm
Brightness	Max 300 cd/m ²
Measuring distance	2.0 ~ 6.0m (0.25m step)
Processor	Amlogic ARM® Cortex®-A5(ARMv7) 1.5Ghz QuadCore Mali™-450 MP2 GPU (OpenGL ES 2.0/1.1)
Memory	SD Memory 8GB 90MB/s
Operating system (OS)	Google's Android
Power supply	DC 12V 5A
Control	Remote controller or Serial communication
Multi-control channel	Max 5 channels
External interface	4x USB, 1x RS-232, 1x Ethernet (10/100)
Standard Accessory	Wall mount bracket
Optional Accessory	Floor stand, Red/Green glasses, Wi-Fi module
Dimension & weight	(W)568 x (H)372 x (D)54mm, Approx. 4.0kg

